


Patient Identifier:	LIVE_COM4	Data Origin:	■ Live Physical Hardware (ESP32/COM4)
Sensor Device:	Physical ESP32 + MPU6050 (COM4)	Sampling Rate:	50.0 Hz (Nyquist: 25.0 Hz)
Assessment Duration:	6.0 seconds	Analysis Window:	3.0 sec (50% sliding overlap)

AI Pattern Classification	Tremor Severity Index	Clinical Severity Grade
HEALTHY Confidence: 99.0% (PD Match: 0.0%)	0.0 / 100	Minimal / Negligible Physiological baseline or voluntary movement. No rhythmic 4-6 Hz tremor detected.

Biomechanical Time-Series & Spectral Isolation

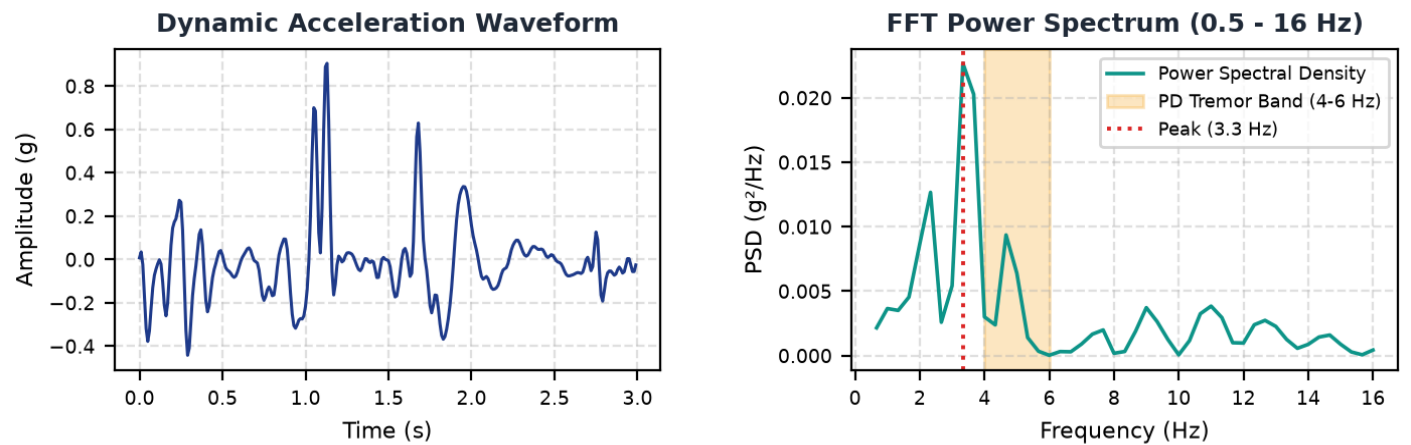

Extracted Digital Biomarkers (Window Metrics)

Biomarker	Observed Value	Reference Band / Diagnostic Context
Dominant Frequency	3.33 Hz	Parkinsonian Resting Tremor: 4.0 - 6.0 Hz
Tremor Band Power Ratio	15.1%	Ratio of power in 4-6 Hz band to broad spectrum (>40% indicates localized tremor)
Dynamic Jerk RMS	7.94 g/s	Time-derivative of acceleration measuring abrupt oscillatory changes
Spectral Shannon Entropy	0.838	Low (<0.45) indicates rhythmic periodicity; High (>0.70) indicates random/voluntary noise

Physiological Interpretation & Explainability

Dominant spectral energy is at 3.3 Hz, reflecting low-frequency voluntary movement or ambient baseline shifts rather than localized tremor. Tremor-band energy is minimal (15.1%), consistent with physiological resting state. Spectral entropy (0.84) reflects broad-spectrum kinetic noise or non-rhythmic movement. Dynamic amplitude RMS (0.185 g) and jerk (7.9 g/s) signify substantial biomechanical displacement and rapid oscillatory acceleration.

NOTICE & MEDICAL DISCLAIMER: NeuroTrack is an experimental screening and longitudinal monitoring aid, NOT a diagnostic or medical treatment device. The metrics, classifications, and severity scores provided herein reflect biomechanical time-series observations and must be interpreted solely by a qualified neurologist or physician in clinical context. Never initiate, alter, or discontinue any medication based on this report.